

Accuracy results summary – Medflow Labs

The following results describe a series of experiments that took place in MedFlow's lab.

For our experiments, we performed 50 sessions of peeing and compared the results that our device calculated to the FlowMaster's results.

The following results show the comparison between to FlowMaster to our device focuses on the 3 main values: Total volume, AVG flow and Max flow.

From the total 50 sessions, 17 of the sessions conducted with urine while the other 33 were conducted using water.

We used to following definitions for accuracy:

Accurate: Results that are close to the FlowMaster with $-5/+5\%$.

Not Accurate: Results that are far from the FlowMaster with more than $-5/+5\%$.

Overall results for the 50 sessions:


For the 50 sessions, 3 sessions found accurate with the Total Volume, 7 sessions found accurate with the AVG flow and 3 sessions found accurate with the Max flow.



Sessions that exceeded by more than 5 percent up: 18 sessions exceeded the Total Volume, 24 sessions exceeded the AVG flow and 20 sessions exceeded the Max flow.

Sessions that exceeded by more than 5 percent down: 29 sessions exceeded the Total Volume, 19 sessions exceeded the AVG flow and 27 sessions exceeded the Max flow.

The following graphs visualities the result:

Color green represents accuracy (close by less than $-/+5\%$), color orange represents sessions exceeded by more than 5 percent up and color red represents sessions exceeded by more than 5 percent down.

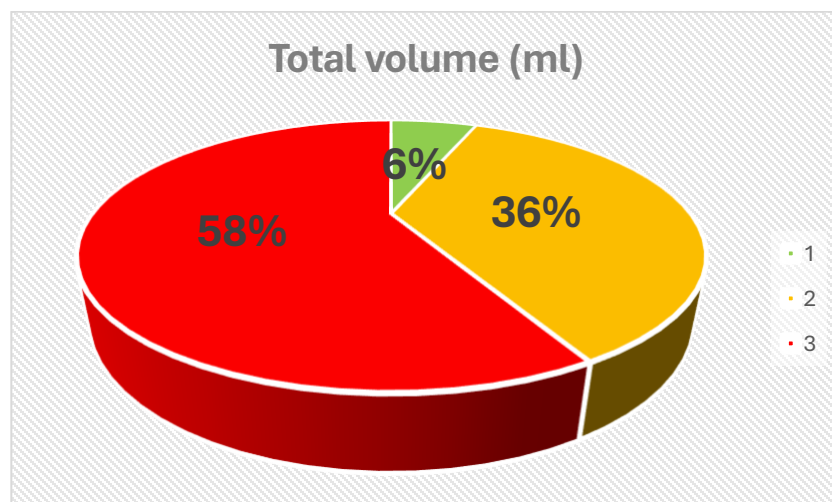

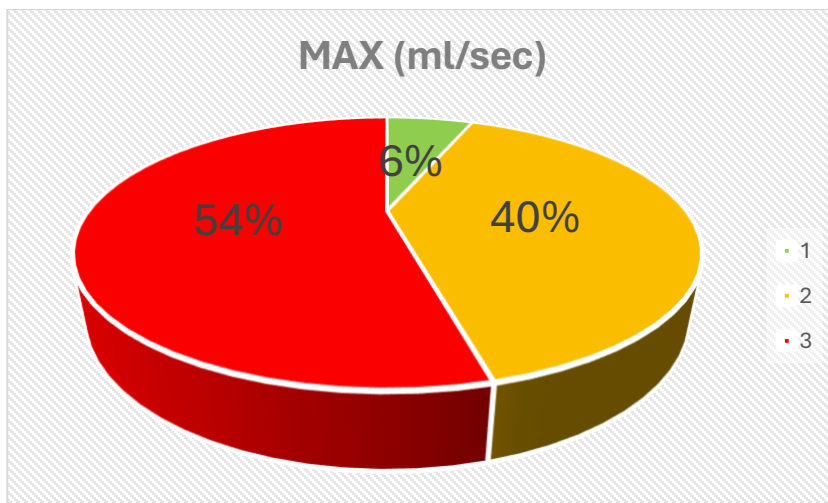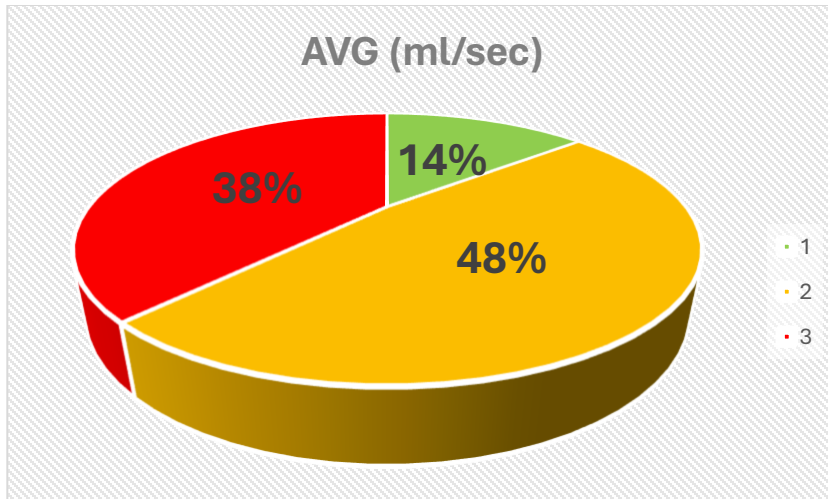

Overall results for the Urine sessions:

From the 50 sessions, 17 have been conducted using urine.

One session has been found accurate with the Total Volume and one session has been found accurate with the Max flow.


Sessions that exceeded by more than 5 percent up: 3 sessions exceeded the Total Volume, 4 sessions exceeded the AVG flow and 6 sessions exceeded the Max flow.

Sessions that exceeded by more than 5 percent down: 13 sessions exceeded the Total Volume, 13 sessions exceeded the AVG flow and 10 sessions exceeded the Max flow.


The following graphs visualities the results:

Color green represents accuracy (close by less than $\pm 5\%$), color orange represents sessions exceeded by more than 5 percent up and color red represents sessions exceeded by more than 5 percent down.


Total volume (ml)

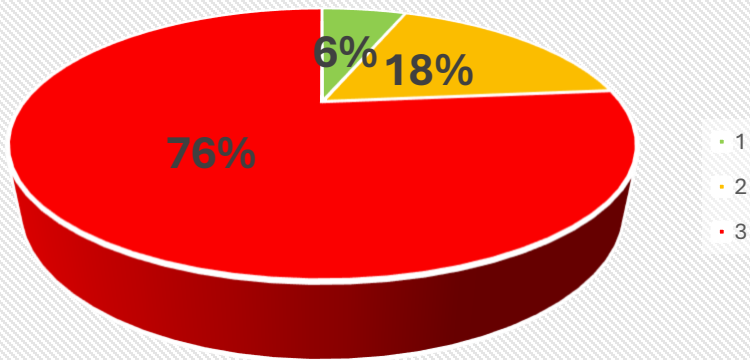

AVG (ml/sec)

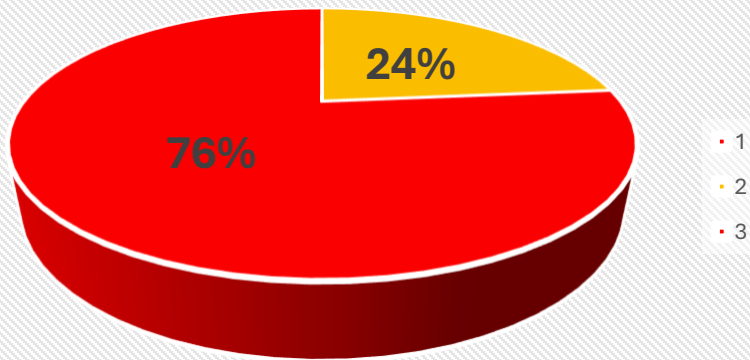

MAX (ml/sec)

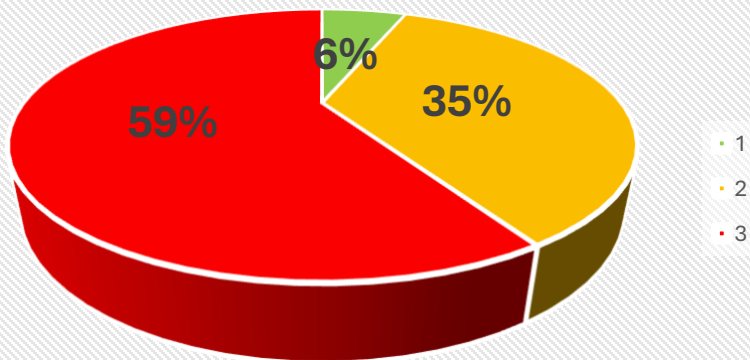

Overall results for the Water sessions:

From the 50 sessions, 33 have been conducted using urine.

2 sessions have been found accurate with the Total Volume, 7 sessions have been found accurate with the AVG flow and 2 sessions have been found accurate with the Max flow

Sessions that exceeded by more than 5 percent up: 15 sessions exceeded the Total Volume, 20 sessions exceeded the AVG flow and 14 sessions exceeded the Max flow.

Sessions that exceeded by more than 5 percent down: 16 sessions exceeded the Total Volume, 6 sessions exceeded the AVG flow and 17 sessions exceeded the Max flow.

The following graphs visualities the results:

Color green represents accuracy (close by less than $\pm 5\%$), color orange represents sessions exceeded by more than 5 percent up and color red represents sessions exceeded by more than 5 percent down.

